

9A279002 Extended Back Section Side Rail

KIT COMPONENTS

Part Number	Description	Qty.
051-0900-00	Extended Side Rail	2
003-1078-00	I / O Instructions.....	1


REQUIRED TOOLS

Phillips Screwdriver

INTRODUCTION

To provide longer rails on the side of the backrest for mounting Midmark approved accessories.

INSTALLATION

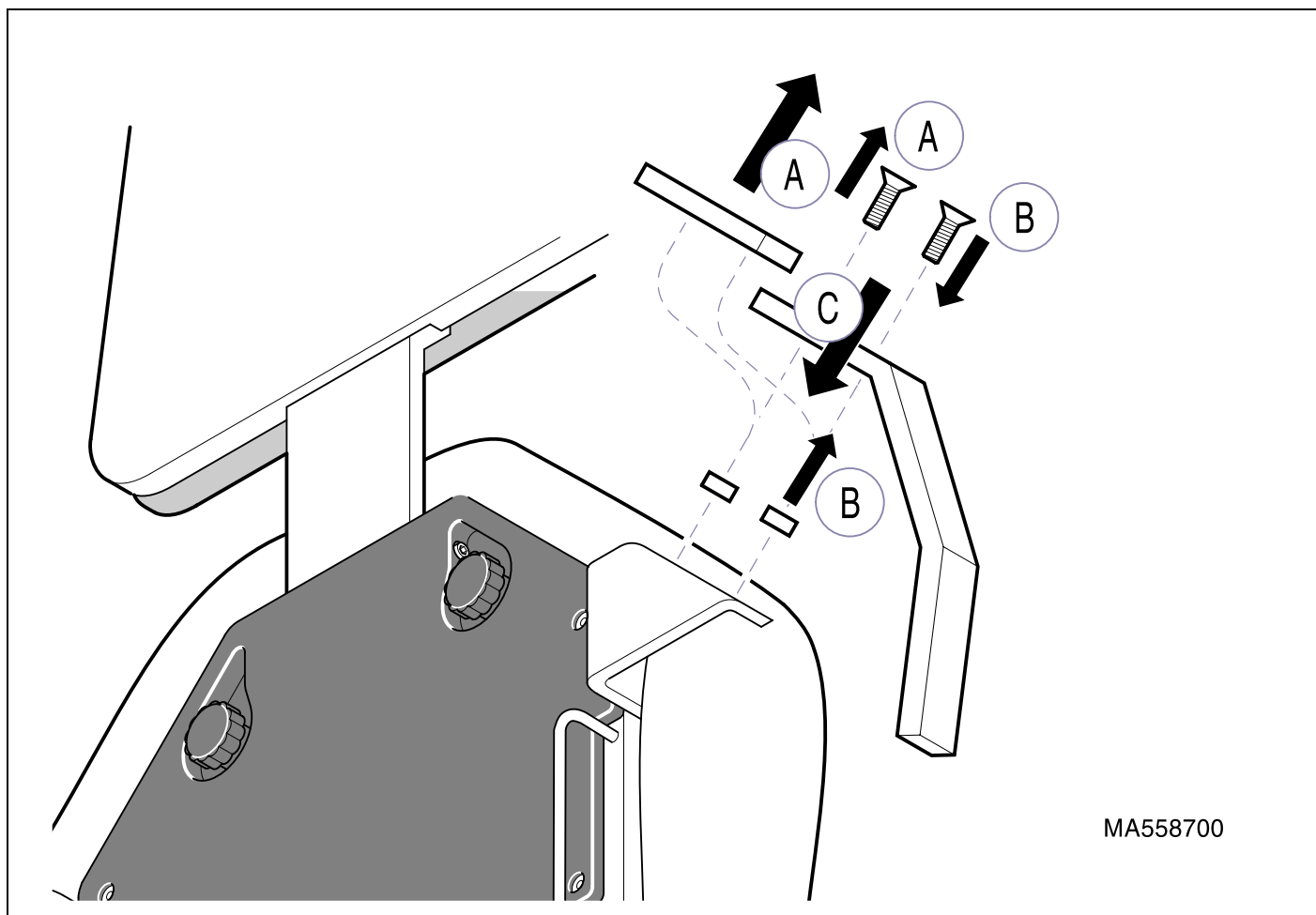

NOTE

Retain the screws and spacers from the existing siderails for use in mounting the new extended siderails.

- A. Remove mounting screws, existing siderails, and spacers.
- B. Insert mounting screws through side rail and slide spacers onto screws.
- C. Secure the side rails and spacers in position on the back rest.